


ottobock.

ParaGolfer

Get up and play

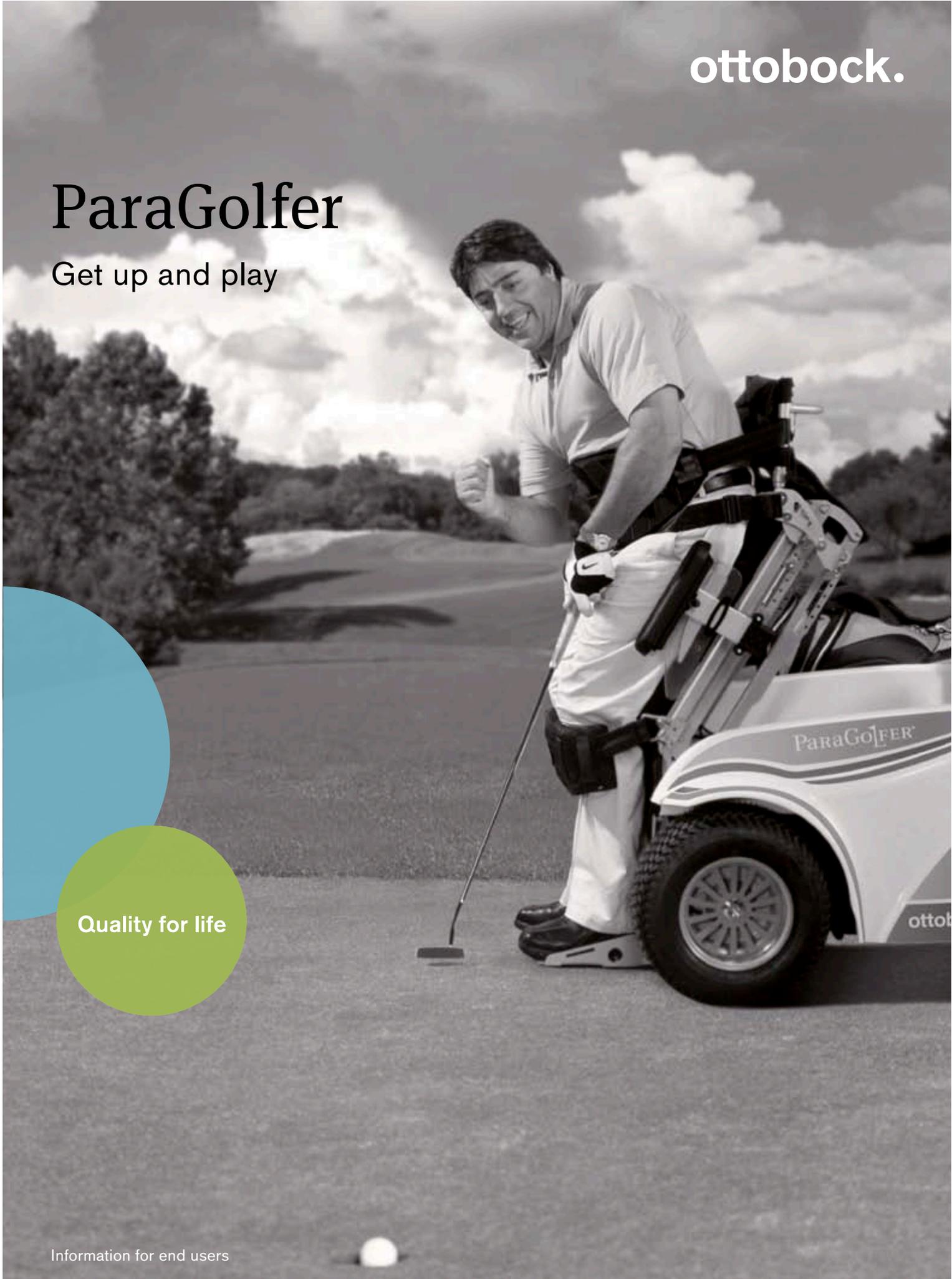

Quality for life

High Quality and Safety



Played in a natural setting, golf is a recreational activity that allows you to spend time with friends and family. But above all, golf is an athletic challenge. With a powerful yet elegant rotation of the upper body, the player brings the club into position, aims and tees off. A unique moment which, thanks to the ParaGolfer, can also be experienced by people with limited mobility. The innovative ParaGolfer by Ottobock offers unparalleled off-road capabilities in combination with a patented stand-up assist and unrivalled stability ...





ParaGolfer

Always in Motion

"Golf means exercise, fun and quality of life. This is especially important for people with limited mobility."

Anthony Netto knows what he is talking about. Even in hospital, the professional golfer suffering from multiple sclerosis recognised the integrative, healing effect of sports activities on the physical and mental well-being of people with physical limitations. Starting to play golf again in spite of his illness was merely a matter of time for Netto.

With the ParaGolfer by Ottobock he has finally reached this goal and is reaping the benefits in all respects: Among other things, he is once again able to hit the ball well over 200 metres. This physical activity speeds up his metabolism, stretches his muscles, reduces spasms and promotes joint mobility. The unique technology of the ParaGolfer brings the player to an upright position, permitting precisely those challenging movements to which Netto owes his 0 handicap.

For years, Anthony Netto has been travelling the world with his ParaGolfer, working with professional golfers and training the German Disabled Golf Team. It is important to him that both able-bodied and disabled players are able to play the game of golf at the same level: His stated goal is to fully reintegrate people with physical limitations into society. You play the ball as it lies, and how good your game is depends solely on the skill of the individual.

Furthermore, since the ParaGolfer puts no greater demands on the carefully manicured lawns than a lawnmower, he and his team are welcome on greens around the world. The manoeuvrable vehicle even allows a player to hit the ball out of a sand trap. "It lets you play a full round of normal golf," Anthony Netto exults.




Anthony Netto

Born in South Africa in 1963, the professional golfer has been wheelchair-bound since the year 2000. With the ParaGolfer, which helps him play golf, Netto supports the German Disabled Golf Team, acts as a worldwide ambassador for handicapped golfers and is the initiator of the "KidSwing" programme.

ParaGolfer

Get up and play

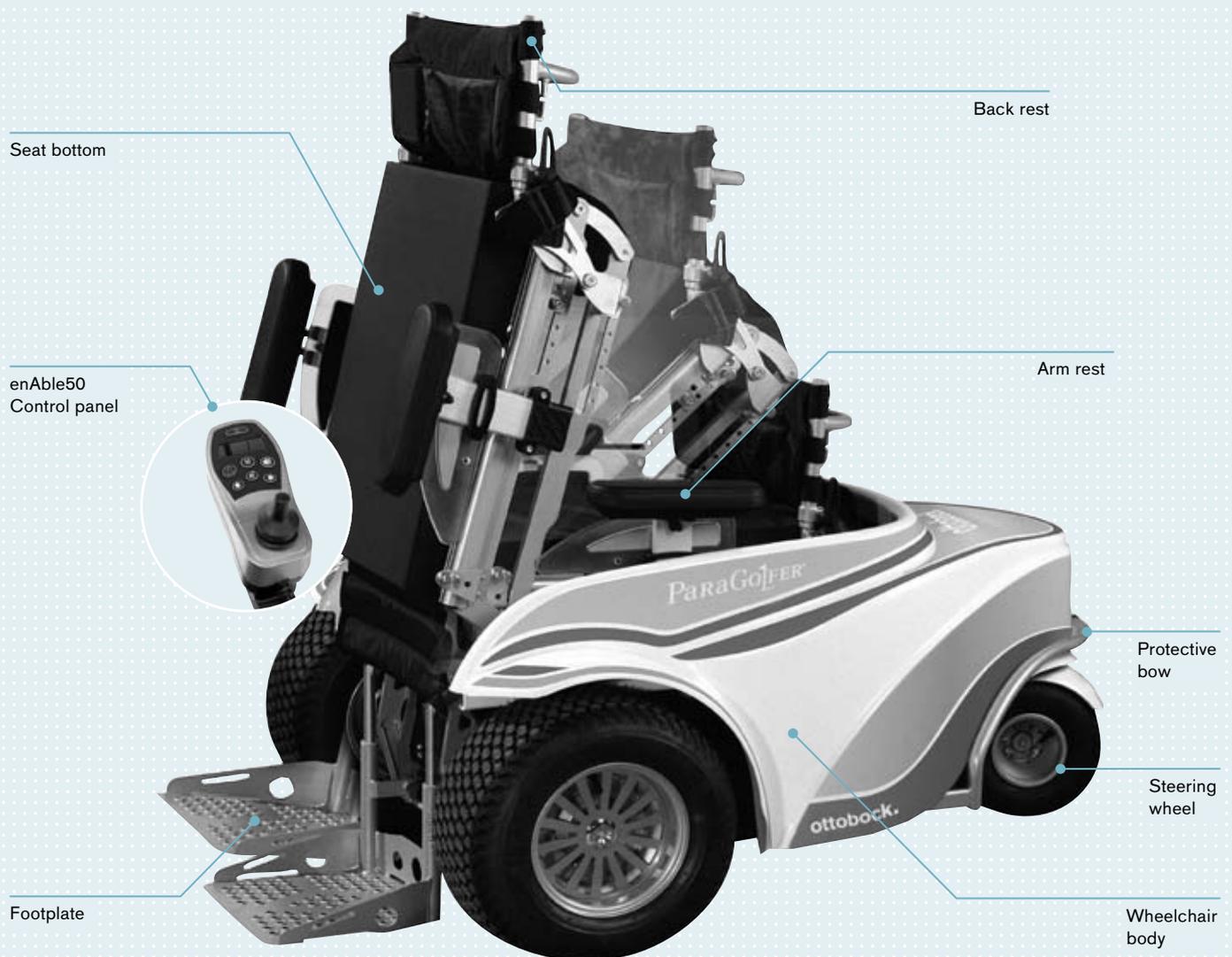

The illustration shows a photomontage with design: Titanium Look

Thanks to its three-wheel design and powerful motors, the ParaGolfer easily conquers slopes of up to 30 % (corresponding to 17°). Neither sand traps nor embankments can stop it. If the maximum values are exceeded, the special off-road design automatically returns the vehicle back to a safe

position, using to the laws of gravity. The player can rely on extremely secure support in any position. This means the freedom to focus on what is really important: Playing the game.



Technical Data	
Seat width	46 cm, can be reduced with lateral pads
Seat depth	36 – 51 cm (2.5 cm increments)
Seat height	52 – 67 cm (continuous)
Seat inclination	0°
Arm rest height	15 cm
Arm rest length	26 cm
Lower leg length	43 – 58 cm (continuous) in combination with the seat height
Back rest height	29 – 35 cm (continuous)
Back angle	+ 5°
Overall width	90 cm
Overall height	From 80 cm depending on back height setting
Overall length	153 cm (driving forward) 149 cm (packing size with folded steering wheel)
Turning radius	190 cm
Empty weight	180 kg
Maximum load capacity	140 kg
Speed	10 km/h
Climbing ability	30 %
Maximum obstacle height	10 cm
Batteries	2x 94 Ah (C5) gel batteries
Range	Approx. 35 km / 18 holes

- The Kids package includes a snap-in footplate for short lower leg lengths, a small seat cushion, a short torso strap and side panels to reduce the seat width.

Otto Bock HealthCare GmbH
Max-Näder-Straße 15, 37115 Duderstadt, Germany
T +49 5527 848-1590 · F +49 5527 848-1676
reha-export@ottobock.de · www.ottobock.com